

Using electron paramagnetic resonance to map N@C60 during high throughput processing

Plant, Simon; Porfyrakis, Kyriakos

DOI:
[10.1039/C4AN00734D](https://doi.org/10.1039/C4AN00734D)

License:
Creative Commons: Attribution (CC BY)

Document Version
Publisher's PDF, also known as Version of record

Citation for published version (Harvard):
Plant, S & Porfyrakis, K 2014, 'Using electron paramagnetic resonance to map N@C₆₀ during high throughput processing', *The Analyst*, vol. 139, no. 18, pp. 4519-4524. <https://doi.org/10.1039/C4AN00734D>

[Link to publication on Research at Birmingham portal](#)

Publisher Rights Statement:
Eligibility for repository : checked 28/01/2015

General rights

Unless a licence is specified above, all rights (including copyright and moral rights) in this document are retained by the authors and/or the copyright holders. The express permission of the copyright holder must be obtained for any use of this material other than for purposes permitted by law.

- Users may freely distribute the URL that is used to identify this publication.
- Users may download and/or print one copy of the publication from the University of Birmingham research portal for the purpose of private study or non-commercial research.
- User may use extracts from the document in line with the concept of 'fair dealing' under the Copyright, Designs and Patents Act 1988 (?)
- Users may not further distribute the material nor use it for the purposes of commercial gain.

Where a licence is displayed above, please note the terms and conditions of the licence govern your use of this document.

When citing, please reference the published version.

Take down policy

While the University of Birmingham exercises care and attention in making items available there are rare occasions when an item has been uploaded in error or has been deemed to be commercially or otherwise sensitive.

If you believe that this is the case for this document, please contact UBIRA@lists.bham.ac.uk providing details and we will remove access to the work immediately and investigate.

CrossMark
click for updatesCite this: *Analyst*, 2014, 139, 4519

Using electron paramagnetic resonance to map N@C₆₀ during high throughput processing

Simon R. Plant†* and Kyriakos Porfyrakis

The endohedral fullerene molecule, N@C₆₀, is a candidate for molecular spin qubits (quantum bits) and spin probes owing to its exceptional electron spin properties. Advancements in the processing of N@C₆₀ are key to obtaining samples of high purity on a reasonable timescale. We investigate enrichment by high throughput processing (flow rate of 18 L h⁻¹ and operating pressure of 1.5–2 MPa) using high performance liquid chromatography (HPLC) as a means of scaling N@C₆₀ production. We use detection by electron paramagnetic resonance (EPR) spectroscopy to map N@C₆₀ during processing, and through the reconstruction of the peak position in the chromatogram, we are able to determine the retention time and relative purity of N@C₆₀ without the need for its isolation. Based on this, we establish a procedure for time-efficient, high throughput processing to isolate N@C₆₀ in high purity.

Received 24th April 2014
Accepted 10th June 2014

DOI: 10.1039/c4an00734d

www.rsc.org/analyst

1 Introduction

High performance liquid chromatography (HPLC) offers a means of separating mixtures of endohedral fullerenes¹ in solution, according to parameters such as mass, structural isomerism and polarizability. The technique is therefore particularly effective in isolating individual endohedral metallofullerenes² owing to pronounced variations in size (mass and structure) of the fullerene cage, as well as variations in the encapsulated elements. It is also possible to isolate endohedral fullerenes such as N@C₆₀ (ref. 3) – which is generated through the direct ion bombardment of C₆₀ – using recycling HPLC,^{4,5} and this is in spite of a yield of 10⁻⁵ to 10⁻⁴ (N@C₆₀/C₆₀) in the starting material and a mass difference of only 2% between the two molecules. N@C₆₀ is paramagnetic, exhibiting narrow EPR (electron paramagnetic resonance) linewidths, and has been selected as a model spin system in which to demonstrate EPR techniques.⁶ Prompted by proposals for a scalable quantum computing architecture,^{7,8} N@C₆₀ has been investigated extensively as a molecular spin qubit (*quantum bit*),^{9–13} utilising the long decoherence time (*e.g.* 0.25 ms at 170 K (ref. 14)) of its electron spin. EPR investigations of N@C₆₀ incorporated into nanostructures,¹⁵ liquid crystals^{16–18} and crystalline matrices^{19,20} have additionally highlighted its potential as a spin probe. N@C₆₀ has also been used for single-molecule transistors,^{21,22} enabling the measurement of spin excitations from electron tunnelling spectra.

In order to avoid the fullerene cluster formation that may influence both the spin lattice relaxation (*T*₁) and phase

decoherence (*T*₂) times, experiments require the preparation of dilute solutions of N@C₆₀ in high purity (>50%). Samples of high purity are also required for schemes that have recently been developed to synthesize N@C₆₀–N@C₆₀ dimers.^{23,24} Such molecules represent isolated pairs of molecular electron spin centres that can permit the investigation of controlled spin–spin coupling and may also enable the demonstration of quantum entanglement of the electron spins. Although there is a demand for super-milligram quantities of N@C₆₀, advancements in production, combined with analysis during the processing of N@C₆₀, are needed. Chromatographic separation with inline EPR detection has been demonstrated before to discriminate between the paramagnetic metallofullerenes Y@C₈₂ and Sc₃@C₈₂.²⁵ HPLC-EPR has the advantage of detecting paramagnetic species amidst empty-cage fullerenes and diamagnetic metallofullerenes, but it has not been widely adopted as a general technique within the field. The drawback is in the selection of parameters for continuous-wave (CW) EPR detection: spectral resolution may be lost at the expense of improved signal-to-noise ratio, a relatively large spectral window may be required to observe the line profile of both species and it is unlikely that a single set of optimal parameters can be used to detect multiple species simultaneously. However, CW-EPR detection is readily applicable to N@C₆₀ as a single set of optimized parameters can be used.

Here, we introduce high-throughput HPLC (flow rate of 18 L h⁻¹) processing of N@C₆₀, which has a capacity of up to 17 times that of previous reports. We use EPR detection to map N@C₆₀ during processing, and through the reconstruction of the peak in the chromatogram, we are able to determine precisely the retention and relative purity of N@C₆₀ as a function of time, without the need for its isolation. This allows us to establish a procedure for time-efficient, high throughput processing in order to isolate N@C₆₀ in high purity.

Department of Materials, University of Oxford, Parks Road, Oxford OX1 3PH, UK.
E-mail: s.r.plant@bham.ac.uk

† Present address: Nanoscale Physics Research Laboratory, School of Physics and Astronomy, University of Birmingham, Birmingham, B15 2TT, UK.

2 Experimental section

C₆₀ powder (99.5+%) was supplied by the MER corporation. The toluene was HPLC grade (99.8%) purchased from Fisher. N@C₆₀ was produced by the ion bombardment of sublimated C₆₀ with ionized nitrogen at a beam energy of 40 eV.²⁶ This method yields a soot consisting of N@C₆₀:C₆₀ in the ratio of typically 1 : 10⁵. The soot was dissolved in toluene and filtered to remove particles greater than 0.2 μm. Batches of fullerene soot containing ¹⁵N@C₆₀ were produced using isotopically enriched nitrogen gas.

The chromatographic separation of N@C₆₀ and C₆₀ was performed using an LC-250HS recycling HPLC apparatus connected to a Cosmosil 15PBB column (100 mm × 500 mm), using injection volumes of 170 mL, which is the maximum loading capacity of the HPLC apparatus. The maximum flow rate of 300 mL min⁻¹ was used throughout, at an operating pressure of 1.5–2 MPa, in order to minimize the retention time of C₆₀ on the column. The UV detector wavelength was fixed at 312 nm during processing. Thereafter, a LC-908W HPLC apparatus connected to a Cosmosil 5PBB column was used for the isolation of high purity N@C₆₀. EPR measurements were made using a continuous-wave (CW) Magnetech Miniscope MS200 spectrometer operating at X-band (~9.5 GHz) and equipped with a rectangular TE₁₀₂ cavity. For comparative studies, the acquisition parameters were identical for all measurements: microwave power of 0.32 mW, modulation amplitude of 0.02 mT, sweep width of 3.50 mT, sweep time of 90 s and 20 scans. This enabled a direct comparison of the relative line intensities of the samples. Optical absorbance was measured using a JASCO V-570 spectrophotometer at a wavelength of 495 nm in ambient conditions.

3 Results and discussion

A comparison of the chromatographic columns used previously for N@C₆₀ enrichment shows that HPLC processing is limited to typical injection volumes of 10 mL and flow rates of 18 mL min⁻¹ at a preparative scale.^{5,27,28} With the aim of increasing throughput, we introduced a new apparatus offering maximum injection volumes of 170 mL and flow rate of 300 mL min⁻¹ coupled to a new chromatographic column with a loading capacity of 1 L. The stationary phase of the chromatographic column consists of 15 μm silica particles functionalized with pentabromobenzyl (PBB) groups. The size of the particles within the stationary phase was selected so as to maintain an operating pressure of 1.5–2 MPa at maximum flow rate. The general procedure is to inject the as-produced mixture of N@C₆₀/C₆₀ into the column, and after the sample elutes, to divide it into two fractions. The first fraction corresponds to the first half of the peak in the chromatogram; the second fraction corresponds to the second half of the peak (as denoted in Fig. 1(a)). Due to the delay in the retention time of N@C₆₀ in the column, the majority of N@C₆₀ within the sample elutes in the second fraction. We can regard fraction 1 as being 'spin-depleted', and fraction 2 as being 'spin-enriched'. Fraction 1 contains (mostly) C₆₀ and is discarded. The remainder of the

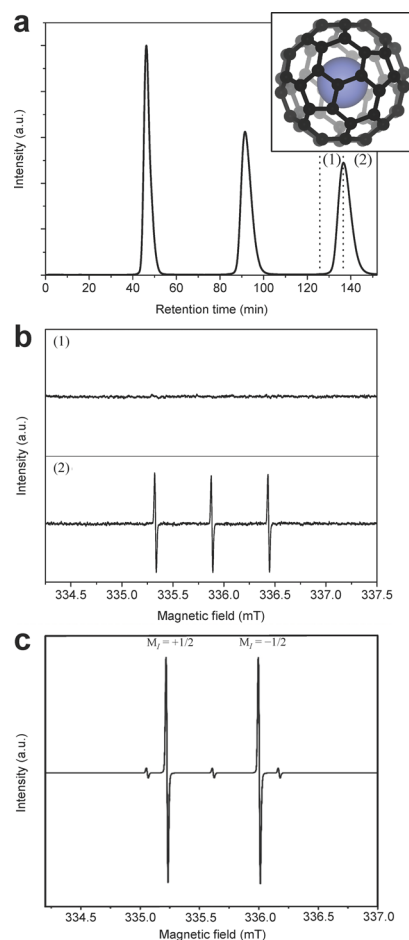


Fig. 1 (a) Chromatogram for a mixture of N@C₆₀ and C₆₀ which passes through the 15PBB column three times and is separated into fractions 1 and 2 on the third cycle. Inset: a model of the N@C₆₀ molecule. (b) Typical CW-EPR spectra at X-band (room temperature) corresponding to fractions 1 and 2, respectively, revealing ¹⁴N@C₆₀ (triplet due to hyperfine splitting) in fraction 2. (c) Example CW-EPR spectrum at X-band showing the hyperfine lines (doublet) for ¹⁵N@C₆₀, dissolved in toluene at room temperature, after several stages of enrichment. The low intensity triplet arises from the presence of a small quantity of ¹⁴N@C₆₀ in the sample.

sample is then injected at time intervals, so that all of the material reaches the same purity before proceeding to stage 2. In stage 2, fraction 2 is injected through the column and the procedure is repeated.

In general, locking the sample in recycling mode during processing by HPLC affords increased separation between the retention times of N@C₆₀ and C₆₀ when eluting from the column. Fig. 1(a) shows a typical chromatogram for an N@C₆₀ and C₆₀ mixture as it passes through the 15PBB column for 3 cycles, whereby the sample is collected during the third cycle. The reasons for collecting the sample during the third cycle were two-fold. (1) It was evident from the chromatogram that the third peak would tail into the fourth due to peak broadening, indicating that the fraction 2 component could become de-enriched after the third cycle. (2) The capacity of the solvent reservoir is 10 L; 7.5 L of toluene were required to collect fractions 1 and 2, shown in Fig. 1(a), at the third cycle. Collecting



the sample on subsequent cycles may well have exhausted the toluene in the reservoir.

This procedure proved effective in containing the N@C_{60} within fraction 2 as evidenced by the EPR spectra shown in Fig. 1(b), which reveals the characteristic triplet of $^{14}\text{N@C}_{60}$. (The more abundant isotope of nitrogen, ^{14}N , was used to produce the N@C_{60} used in these initial studies.) At this flow rate, the retention time of C_{60} is *circa* 45 minutes. Consequently, one injection lasting 3 cycles takes ~ 152 minutes (as shown in Fig. 1). However, using the LC-250HS recycling HPLC apparatus, we are able to inject up to 170 mL at a time. Each sample has to pass through 6 stages of HPLC processing before it is pure enough to move to the next phase of the process.

As $^{15}\text{N@C}_{60}$ was used in further work to demonstrate mapping in the chromatogram, Fig. 1(c) shows an example of an X-band CW-EPR spectrum of dilute $^{15}\text{N@C}_{60}$ (dissolved in toluene and at room temperature) following several stages of enrichment. The $^{15}\text{N@C}_{60}$ spin system consists of an $S = 3/2$ electron spin coupled to an $I = 1/2$ nuclear spin. The spectrum reveals the hyperfine lines (intense doublet) corresponding to the $M_I = +1/2$ and $M_I = -1/2$ nuclear spin states. The hyperfine lines of $^{14}\text{N@C}_{60}$ (low intensity triplet), corresponding to the $M_I = \pm 1, 0$ nuclear spin states of ^{14}N , are visible above the baseline. The presence of $^{14}\text{N@C}_{60}$ is due to the ^{14}N content in the nitrogen gas when $^{15}\text{N@C}_{60}$ is first synthesized.

3.1 Mapping N@C_{60} in the chromatogram

Due to the close retention times of N@C_{60} and C_{60} , it is not possible to resolve the peak corresponding N@C_{60} in the chromatogram directly during early processing. However, by collecting a larger number of sub-fractions and measuring the EPR spectra of a standard aliquot from each fraction, the N@C_{60} peak can be reconstructed and its position mapped. To this end, 170 mL of an as-produced $^{15}\text{N@C}_{60}/\text{C}_{60}$ mixture was injected into the 15PBB column, and after one cycle, the fraction corresponding to the second half of the peak was split into sub-fractions, where each sub-fraction was collected over a period of one minute (see Fig. 2). In a second experiment, we injected another 170 mL of the mixture into the column, and locked the sample in recycling mode, so that the sample passed through the column for a total of 3 cycles (also shown in Fig. 2). On the third cycle, we repeated the procedure for the first experiment, collecting 10 sub-fractions, where each sub-fraction was collected over a period of 1.5 minutes (with the exception of sub-fraction iv, which was collected over a period of 1 minute).

An aliquot from each sub-fraction was taken and concentrated ten times. 0.15 mL of the solution was placed in a quartz tube, in order to carry out EPR measurements. As for Fig. 1(c), the most intense lines in the EPR spectra arise from $^{15}\text{N@C}_{60}$, with the 3 lines corresponding to $^{14}\text{N@C}_{60}$ observable just above the level of the baseline noise. These weak lines arising from $^{14}\text{N@C}_{60}$ were disregarded for the purpose of this experiment; measurements were made on the lines corresponding to $^{15}\text{N@C}_{60}$ only. The line intensity was determined by measuring the peak-to-trough amplitude of each of the two hyperfine lines in the spectrum and taking the average. As the linewidths are

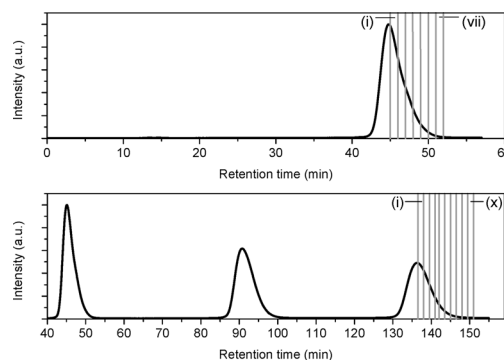


Fig. 2 Chromatograms showing the division of the 'spin-rich' fraction of a $\text{N@C}_{60}/\text{C}_{60}$ mixture into sub-fractions after one cycle (top) and three cycles (bottom). The first and last sub-fractions are marked.

invariant between spectra, measuring the amplitude of the hyperfine lines offers a convenient method for quantitative analysis,²⁹ because under this condition, the signal amplitude is proportional to the number of spins measured. Fig. 3(a) shows a plot of the line intensity as a function of time elapsed in the HPLC chromatogram for the 1st and 3rd cycles through the 15PBB column. This effectively maps the position of the peak corresponding to N@C_{60} in the chromatogram at cycles 1 and 3. The bars in the plot denote the sub-fractions corresponding to the chromatograms shown in Fig. 2. The data have been fitted with modified (skewed) Gaussian distributions, which are otherwise used for fitting chromatograms. It is evident from Fig. 3(a), in comparing the maxima of the reconstructed peaks, that the peak of N@C_{60} has shifted to the right by a time of 3 min 18 s after 3 cycles, while the intensity falls by approximately two-thirds. These data also indicate elution times of 1 min 18 s and 4 min 36 s (relative to the C_{60} peak in the chromatogram) for cycles 1 and 3, respectively. This compares favourably with the data for other chromatographic columns used for N@C_{60} processing on a preparative scale. Taking the ratio of the retention time of N@C_{60} to that of C_{60} provides a simple measure of the efficiency of the column. Based on reports for the 5PYE and 5PBB columns,⁵ these ratios are approximately 1.01 and 1.03, respectively. The ratio for the 15PBB column in the present work is 1.03 (based on analysis for both cycles 1 and 3), indicating that the column exhibits comparable efficiency in separating N@C_{60} from C_{60} , based on peak maxima.

Fig. 3(a) also confirms that the majority of N@C_{60} is contained in the second fraction. However for both cycles 1 and 3, a proportion of the N@C_{60} is retained within the first fraction. One motivation for this study was to assess the efficacy of collecting fraction 2 on the third cycle, rather than on the first. Whilst the separation in the elution times for N@C_{60} and C_{60} is enhanced at the third cycle (as measured by the peak maxima), the distribution changes due to peak broadening, and consequently its intensity is diminished. Significantly, the areas under the curves, proportional to the number of electron spins (and therefore to the quantity of N@C_{60}) are comparable for cycles 1 and 3. Therefore, if one wishes to obtain purified N@C_{60} time-efficiently whilst keeping losses to a minimum, then it is recommended that one cycle through the 15PBB column is



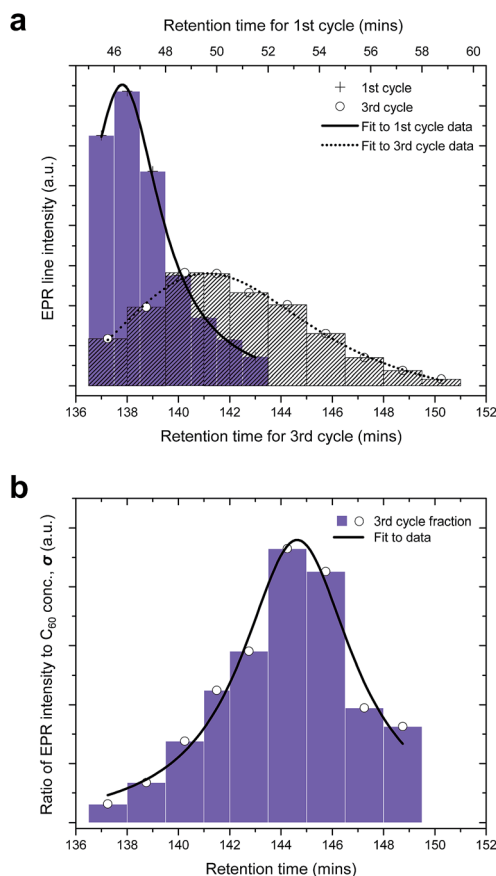


Fig. 3 Mapping the position of the N@C₆₀ peak in the chromatogram. (a) Plot showing EPR line intensity as a function of time elapsed in the HPLC chromatogram after cycles 1 and 3. The bars indicate the fractions collected corresponding to Fig. 2. (b) The evolution of the ratio of the EPR signal intensity and the concentration of C₆₀ (σ) as a function of time elapsed in the HPLC chromatogram.

sufficient in future work. Alternatively, smaller quantities of higher purity N@C₆₀ can be collected during the third cycle, as discussed below.

As a measure of sample purity, we introduce a figure of merit, σ , which is the relative spin concentration (spin centres with respect to C₆₀). We determine σ by taking the ratio of the EPR line intensity to the concentration of C₆₀, given that the line intensity is proportional to the number of spins and molar concentration of the sample is proportional to the number of C₆₀ molecules. In order to calculate the concentration of C₆₀ in each of the sub-fractions, several samples of C₆₀ in HPLC grade toluene were prepared at known concentrations and the absorbance of each of these samples was measured. The molar absorption coefficient, ϵ , was extracted by plotting a calibration curve of absorbance against concentration for the aforementioned samples. The absorbance (at 495 nm) of each of the sub-fractions was measured, and the concentrations were calculated using Beer–Lambert's Law, $A = \epsilon cl$, where A is the absorbance, ϵ is the molar absorption coefficient (as stated above), c is the molar concentration and l is the pathlength. Fig. 3(b) shows a plot of σ as a function of the time elapsed in the HPLC chromatogram for cycle 3. The best fit to the data was achieved using

a Lorentzian distribution. Fig. 3(b) reveals that the highest purity of N@C₆₀ can be obtained from sub-fractions vi and vii, by collecting between *circa* 143 and 147 minutes.

3.2 Isolation of ¹⁵N@C₆₀ in high purity

Two batches of ¹⁵N@C₆₀/C₆₀ mixtures were prepared using the 15PBB column. At the end of stage 6, the total volume of solution had reduced to <20 mL, which is <12% of the total injection capacity permitted by the LC-250HS recycling HPLC apparatus. Therefore, after stage 6, it is no longer advantageous to use the 15PBB for enrichment. It was appropriate to move to the LC-908W HPLC apparatus connected to a 5PBB column, for which the injection volume capacity is 10 mL and the retention time for C₆₀ is ~11 minutes for a flow rate of 18 mL min⁻¹ at ambient temperature. Fig. 4 shows a typical chromatogram for a ¹⁵N@C₆₀/C₆₀ mixture passing through the 5PBB column (eluent: toluene, flow rate: 18 mL min⁻¹, injection volume: 10 mL, ambient temperature). Firstly, minor impurities appear with retention times *circa* 4 and 9 minutes, which arise following concentration of the sample using a rotary evaporator, and these are discarded. The creation of these impurities may be limited by reducing the bath temperature of the rotary evaporator to <30 °C, which also reduces rate of evaporation thereby increasing the time required to concentrate samples. Recycling of the ¹⁵N@C₆₀/C₆₀ mixture through the 15PBB column is used in order to monitor the emergence of other impurities – an example being N₂@C₆₀ (ref. 28) – as marked by the arrow in the inset of Fig. 4. Such impurities were not removed, but were collected along with the second fraction when the sample was divided into two fractions (corresponding to the first and second halves of the peak) during the sixth cycle. N₂@C₆₀ was removed from the ¹⁵N@C₆₀/C₆₀ mixture using recycling HPLC at a later stage.

The chromatogram from the final stage of processing is shown in Fig. 5(a). Recycling HPLC was used to separate C₆₀ from ¹⁵N@C₆₀ (eluent: toluene; flow rate: 18 mL min⁻¹; column: 5PBB). The single peak in the chromatogram begins to split into two, corresponding to C₆₀ and ¹⁵N@C₆₀, during the second

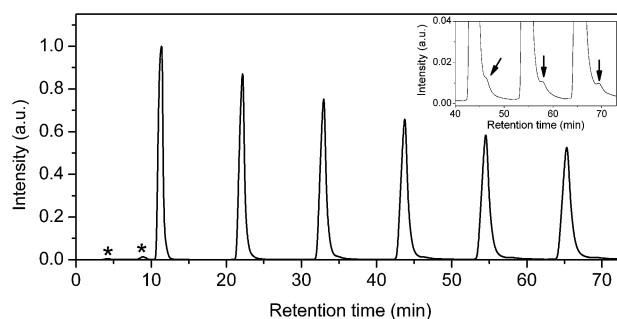


Fig. 4 A typical HPLC chromatogram for a ¹⁵N@C₆₀/C₆₀ mixture as it passes through the 5PBB column. (*) denotes minor impurities that appear following concentration of the sample using rotary evaporation. Inset: amplification ($\times 25$) of the signal intensity for the aforementioned chromatogram peak. The arrow indicates the emergence of an impurity from the main C₆₀ peak during cycles 4–6.



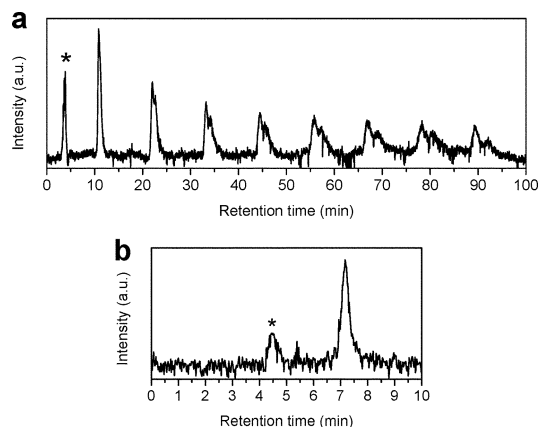


Fig. 5 Isolation of a high purity $^{15}\text{N@C}_{60}$ sample using HPLC. (a) Separation of C_{60} and $^{15}\text{N@C}_{60}$, which appear as two peaks in the chromatogram during the final stage of processing by recycling HPLC (eluent: toluene; flow rate: 18 mL min^{-1} ; column: 5PBB). $^{15}\text{N@C}_{60}$ corresponds to the second peak of the doublet. (b) Final pass through the Buckyprep-M column to permit the quantification of the sample (eluent: toluene; flow rate: 18 mL min^{-1}).

cycle. At the eighth cycle, two fractions were collected corresponding to these peaks. Due to the way the peaks corresponding to C_{60} and $^{15}\text{N@C}_{60}$ overlap at the eighth cycle, this provided an $^{15}\text{N@C}_{60}$ sample of >50% purity. This high purity $^{15}\text{N@C}_{60}$ sample was then concentrated and injected through the Buckyprep-M column (eluent: toluene; flow rate: 18 mL min^{-1}) in order to quantify the sample based on the integrated peak area (see Fig. 5(b)). Yields achieved by following the procedure set out here will be comparable with other procedures for the HPLC processing of N@C_{60} .

4 Conclusions

We have introduced high throughput HPLC processing for the purification (enrichment) of N@C_{60} , using flow rates of 18 L hour^{-1} and pressures of up to 2 MPa. The use of a high capacity HPLC apparatus permits injection volumes of 170 mL, which is up to $17\times$ the injection volumes used in previous work. Further modifications of the HPLC apparatus may enable us to take advantage of the maximum loading capacity of the 15PBB column by injecting 1 L of sample solution per injection.

Additionally, we have used electron paramagnetic resonance spectroscopy as a means to map N@C_{60} during high throughput processing, in order to determine the retention times and relative purity without the need for its isolation. Mapping the peak in the chromatogram has enabled us to demonstrate that the 15PBB column itself has a comparable efficiency to the 5PYE and 5PBB columns used previously, and it has allowed us to develop a procedure for the time-efficient, high throughput enrichment of N@C_{60} .

The goal of the present work was to obtain the largest absolute quantities of N@C_{60} at the highest possible purity on a reasonable timescale. Therefore, we conclude the most time-efficient means of purifying N@C_{60} whilst minimizing losses is to

use the high capacity apparatus to pass a sample containing N@C_{60} through the 15PBB column for one cycle, and then to collect the fraction corresponding to the second half of the peak in the chromatogram. The process is then repeated 6 times until the same sample has been enriched sufficiently to move to a smaller column, such as the 5PBB, whereby the sample is injected through the column in recycling mode, ultimately permitting the isolation of high purity N@C_{60} . Alternatively, we have shown that, in order to collect smaller quantities of high purity N@C_{60} efficiently (*without the constraint to minimize losses*), the following route is possible. The sample containing N@C_{60} is injected into the high capacity HPLC system, locked in recycling mode for 2 cycles, and then at the 3rd cycle, the N@C_{60} -rich sample is collected between *circa* 143 and 147 minutes. The collected sample is concentrated and the procedure is repeated until the required purity has been reached.

The method of scaled processing demonstrated here improves the feasibility of utilising N@C_{60} for applications, such as for spin probes or molecular qubits.

Acknowledgements

We thank Andrew Briggs for supporting this work, and Hippolyte Queudet for his assistance. S. R. P. acknowledges support from EPSRC. K. P. is also supported by EPSRC (EP/K030108/1).

References

- 1 A. A. Popov, S. Yang and L. Dunsch, *Chem. Rev.*, 2013, **113**, 5989–6113.
- 2 H. Shinohara, *Rep. Prog. Phys.*, 2000, **63**, 843–892.
- 3 T. Almeida Murphy, Th. Pawlik, A. Weidinger, M. Höhne, R. Alcalá and J.-M. Spaeth, *Phys. Rev. Lett.*, 1996, **77**, 1075–1078.
- 4 P. Jakes, K.-P. Dinse, C. Meyer, W. Harneit and A. Weidinger, *Phys. Chem. Chem. Phys.*, 2003, **5**, 4080–4083.
- 5 M. Kanai, K. Porfyrakis, G. A. D. Briggs and T. J. S. Dennis, *Chem. Commun.*, 2004, 210–211.
- 6 J. J. L. Morton, A. M. Tyryshkin, A. Ardavan, K. Porfyrakis, S. A. Lyon and G. A. D. Briggs, *J. Chem. Phys.*, 2005, **122**, 174504.
- 7 W. Harneit, *Phys. Rev. A: At., Mol., Opt. Phys.*, 2002, **65**, 032322.
- 8 S. C. Benjamin, A. Ardavan, G. A. D. Briggs, D. A. Britz, D. Gunlycke, J. Jefferson, M. A. G. Jones, D. F. Leigh, B. W. Lovett, A. N. Khlobystov, S. A. Lyon, J. J. L. Morton, K. Porfyrakis, M. R. Sambrook and A. M. Tyryshkin, *J. Phys.: Condens. Matter*, 2006, **18**, S867–S883.
- 9 C. Meyer, W. Harneit, B. Naydenov, K. Lips and A. Weidinger, *Appl. Magn. Reson.*, 2004, **27**, 123–132.
- 10 J. J. L. Morton, A. M. Tyryshkin, A. Ardavan, K. Porfyrakis, S. A. Lyon and G. A. D. Briggs, *Phys. Rev. Lett.*, 2005, **95**, 200501.
- 11 J. J. L. Morton, A. M. Tyryshkin, A. Ardavan, S. C. Benjamin, K. Porfyrakis, S. A. Lyon and G. A. D. Briggs, *Nat. Phys.*, 2006, **2**, 40–43.



- 12 G. W. Morley, J. van Tol, A. Ardavan, K. Porfyakis, J. Zhang and G. A. D. Briggs, *Phys. Rev. Lett.*, 2007, **98**, 220501.
- 13 R. M. Brown, A. M. Tyryshkin, K. Porfyakis, E. M. Gauger, B. W. Lovett, A. Ardavan, S. A. Lyon, G. A. D. Briggs and J. J. L. Morton, *Phys. Rev. Lett.*, 2011, **106**, 110504.
- 14 J. J. L. Morton, A. M. Tyryshkin, A. Ardavan, K. Porfyakis, S. A. Lyon and G. A. D. Briggs, *J. Chem. Phys.*, 2006, **124**, 014508.
- 15 T. Wakahara, T. Kato, K. Miyazawa and W. Harneit, *Carbon*, 2012, **50**, 1709–1712.
- 16 P. Jakes, N. Weiden, R.-A. Eichel, A. Gembus, K.-P. Dinse, C. Meyer, W. Harneit and A. Weidinger, *J. Magn. Reson.*, 2002, **156**, 303–308.
- 17 C. Meyer, W. Harneit, K. Lips, A. Weidinger, P. Jakes and K.-P. Dinse, *Phys. Rev. A: At., Mol., Opt. Phys.*, 2002, **65**, 061201.
- 18 G. Liu, M. C. Gimenez-Lopez, M. Jevric, A. N. Khlobystov, G. A. D. Briggs and K. Porfyakis, *J. Phys. Chem. B*, 2013, **117**, 5925–5931.
- 19 B. Naydenov, C. Spudat, W. Harneit, H. I. Süss, J. Hulliger, J. Nuss and M. Jansen, *Chem. Phys. Lett.*, 2006, **424**, 327–332.
- 20 J. Yang, P. Feng, A. Sygula, W. Harneit, J.-H. Su and J. Du, *Phys. Lett. A*, 2012, **376**, 1748–1751.
- 21 N. Roch, R. Vincent, F. Elste, W. Harneit, W. Wernsdorfer, C. Timm and F. Balestro, *Phys. Rev. B: Condens. Matter Mater. Phys.*, 2011, **83**, 081407(R).
- 22 J. E. Grose, E. S. Tam, C. Timm, M. Scheloske, B. Ulgut, J. J. Parks, H. D. Abruña, W. Harneit and D. C. Ralph, *Nat. Mater.*, 2008, **7**, 884–889.
- 23 B. J. Farrington, M. Jevric, G. A. Rance, A. Ardavan, A. N. Khlobystov, G. A. D. Briggs and K. Porfyakis, *Angew. Chem., Int. Ed.*, 2012, **51**, 3587–3590.
- 24 S. R. Plant, M. Jevric, J. J. L. Morton, A. Ardavan, A. N. Khlobystov, G. A. D. Briggs and K. Porfyakis, *Chem. Sci.*, 2013, **4**, 2971–2975.
- 25 S. Stevenson, H. C. Dorn, P. Burbank, K. Harich, Z. Sun, C. H. Kiang, J. R. Salem, M. S. DeVries, P. H. M. van Loosdrecht, R. D. Johnson, C. S. Yannoni and D. S. Bethune, *Anal. Chem.*, 1994, **66**, 2680–2685.
- 26 K. Lips, M. Waiblinger, B. Pietzak and A. Weidinger, *Phys. Status Solidi A*, 2000, **177**, 81–91.
- 27 A. Weidinger, M. Waiblinger, B. Pietzak and T. Almeida Murphy, *Appl. Phys. A: Mater. Sci. Process.*, 1998, **66**, 287–292.
- 28 T. Suetsuna, N. Dragoe, W. Harneit, A. Weidinger, H. Shimotani, S. Ito, H. Takagi and K. Kitazawa, *Chem.–Eur. J.*, 2002, **8**, 5079–5083.
- 29 G. R. Eaton, S. S. Eaton, D. P. Barr and R. T. Weber, *Quantitative EPR*, Springer, Vienna, 2010.

