

Generation of human induced pluripotent stem cells in defined, feeder-free conditions

Kyrkou, Athena; Stellas, Dimitris; Syrrou, Marika; Klinakis, Apostolos; Fotsis, Theodore; Murphy, Carol

DOI:

[10.1016/j.scr.2016.05.006](https://doi.org/10.1016/j.scr.2016.05.006)

License:

Creative Commons: Attribution-NonCommercial-NoDerivs (CC BY-NC-ND)

Document Version

Publisher's PDF, also known as Version of record

Citation for published version (Harvard):

Kyrkou, A, Stellas, D, Syrrou, M, Klinakis, A, Fotsis, T & Murphy, C 2016, 'Generation of human induced pluripotent stem cells in defined, feeder-free conditions', *Stem Cell Research*.
<https://doi.org/10.1016/j.scr.2016.05.006>

[Link to publication on Research at Birmingham portal](#)

General rights

Unless a licence is specified above, all rights (including copyright and moral rights) in this document are retained by the authors and/or the copyright holders. The express permission of the copyright holder must be obtained for any use of this material other than for purposes permitted by law.

- Users may freely distribute the URL that is used to identify this publication.
- Users may download and/or print one copy of the publication from the University of Birmingham research portal for the purpose of private study or non-commercial research.
- User may use extracts from the document in line with the concept of 'fair dealing' under the Copyright, Designs and Patents Act 1988 (?)
- Users may not further distribute the material nor use it for the purposes of commercial gain.

Where a licence is displayed above, please note the terms and conditions of the licence govern your use of this document.

When citing, please reference the published version.

Take down policy


While the University of Birmingham exercises care and attention in making items available there are rare occasions when an item has been uploaded in error or has been deemed to be commercially or otherwise sensitive.

If you believe that this is the case for this document, please contact UBIRA@lists.bham.ac.uk providing details and we will remove access to the work immediately and investigate.



Contents lists available at ScienceDirect

Stem Cell Research

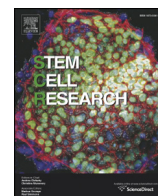journal homepage: www.elsevier.com/locate/scr

Lab resource

Generation of human induced pluripotent stem cells in defined, feeder-free conditions

Athena Kyrkou^a, Dimitris Stellas^b, Marika Syrrou^c, Apostolos Klinakis^b, Theodore Fotsis^{a,d,e}, Carol Murphy^{a,e,*}^a Department of Biomedical Research, Institute of Molecular Biology and Biotechnology, Foundation of Research and Technology-Hellas, University Campus, 45110 Ioannina, Greece^b Centre for Basic Research, Biomedical Research Foundation of the Academy of Athens, 11527 Athens, Greece^c Laboratory of Biology, Medical School, University of Ioannina, 45110 Ioannina, Greece^d Laboratory of Biological Chemistry, Medical School, University of Ioannina, 45110 Ioannina, Greece^e School of Biosciences, College of Life and Environmental Sciences, University of Birmingham, Edgbaston, Birmingham, B15 2TT, UK

ARTICLE INFO

Available online xxxx

ABSTRACT

Herein, we describe a modified protocol for the generation of human induced pluripotent stem cells (hiPS) and expansion under defined, serum free and feeder free conditions. These cells exhibit a high level of plasticity towards various differentiation pathways both in vitro and in vivo. Ultimately, hiPS-derived lines achieved high standards of three dimensional differentiations on biomaterial scaffolds and promoted in vivo regeneration of complex organs, such as Anterior Cruciate Ligament (in swine ACL-rupture models) and other tissues as well.

© 2016 Published by Elsevier B.V. This is an open access article under the CC BY-NC-ND license (<http://creativecommons.org/licenses/by-nc-nd/4.0/>).

Resource table.

Name of stem cell construct	FORTHi001A (hiPS line 1) FORTHi001B (hiPS line 2) FORTHi001C (hiPS line 3)
Institution	Department of Biomedical Research- Institute of Molecular Biology & Biotechnology, Foundation of Research and Technology-Hellas (BRI-FORTH), University Campus, 45110 Ioannina, Greece
Person who created resource	Athena Kyrkou
Contact person and email	Carol Murphy, c.e.murphy@bham.ac.uk
Date archived/stock date	30 June 2014
Origin	Human skin fibroblasts
Type of resource	Human induced pluripotent cells (hiPS)
Sub-type	Cell line
Key transcription factors	Human Oct3/4, Sox2, L-Myc, Klf4, Lin28, sh-p53
Authentication	Identity and purity of cell line confirmed (Fig. 1)
Link to related literature (direct URL links and full references)	Kouroupis et al., 2016, Stem Cell Research-submitted
Information in public databases	Microarray data available in NCBI Gene Expression Omnibus http://www.ncbi.nlm.nih.gov/geo/GEO#GSE58932

DOI of original article: <http://dx.doi.org/10.1016/j.scr.2016.04.016>.

* Corresponding author at: University Campus, Department of Biomedical Research, Institute of Molecular Biology & Biotechnology, Foundation of Research and Technology-Hellas, 45110 Ioannina, Greece.

E-mail address: kouroupis_d@yahoo.gr (T. Fotsis).

1. Resource details, materials and methods

Human skin fibroblasts (ATCC #CRL 2429) were plated at a density of 100,000 cells per 35 mm cell culture dish. The reprogramming strategy was based on a modification of the Y4 episomal approach (Okita et al.). Plasmids for the reprogramming pCXLE-hOCT3/4 -shp53pCXLE -hULpCXLE were purchased from Addgene, Middlesex, UK. In our modified protocol, subconfluent fibroblasts were transfected with the episomes using lipofectamine 2000 (LifeSciences, #11668, Invitrogen, Carlsbad, CA, USA), according to the manufacturers' protocol. 48 h post-transfection, cells were refreshed with fibroblast medium (IMDM #30–2005, supplemented with 10% fetal bovine serum/FBS #16000–044, Gibco, Gaithersburg, USA). On day 4 cells were fed with serum containing embryonic stem cell medium (Knockout DMEM/F12 #12660, 15% FBS #16000–044, 1% Non-Essential Amino Acids/NEAA #11140, 1% Glutamax #35050, 1 × 2-mercaptoethanol #21985023 all from Gibco, supplemented with 40 ng/ml bfgf #100–18B, Peprotech, Rocky Hill, USA) which was refreshed daily, for the next two days. On day 6, this cell population was trypsinized and reseeded 1:7 on irradiated fibroblasts (#GSC-6301G, GlobalStem, Oxfordshire, UK) to boost reprogramming. On day 8 cells were grown in serum free embryonic stem cell medium (Knockout DMEM/F12 #12660, 15% KnockOut Serum Replacement #10828 -Gibco, 1% NEAA, 1% Glutamax, supplemented with 1 × 2-mercaptoethanol and 20 ng/ml bfgf) and the medium was changed daily. From day 12 to day 20 cells were fed with conditioned serum-free embryonic stem cell medium from feeders, to sustain pre-formed human induced pluripotent (pre-hiPS) colonies' growth. On day 20, pre-hiPS

<http://dx.doi.org/10.1016/j.scr.2016.05.006>1873–5061/© 2016 Published by Elsevier B.V. This is an open access article under the CC BY-NC-ND license (<http://creativecommons.org/licenses/by-nc-nd/4.0/>).

Please cite this article as: Kyrkou, A., et al., Generation of human induced pluripotent stem cells in defined, feeder-free conditions, Stem Cell Res. (2016), <http://dx.doi.org/10.1016/j.scr.2016.05.006>

were manually picked using 1 mg/ml dispase (#17105, Gibco) and grown on matrigel (dilution 1:100, # 354277, BD Biosciences, Oxford, UK) with mTeSR medium (#05850, Stem Cell Technologies, Vancouver, Canada). Thereafter, hiPS colonies were manually dissected “en block” and expanded in feeder free and serum free conditions (Fig. 1A).

2. Verification and authentication

2.1. Phenotypic characterization of hiPS

Successful reprogramming of initial fibroblasts towards induced pluripotency was evaluated by immunostaining with key pluripotency

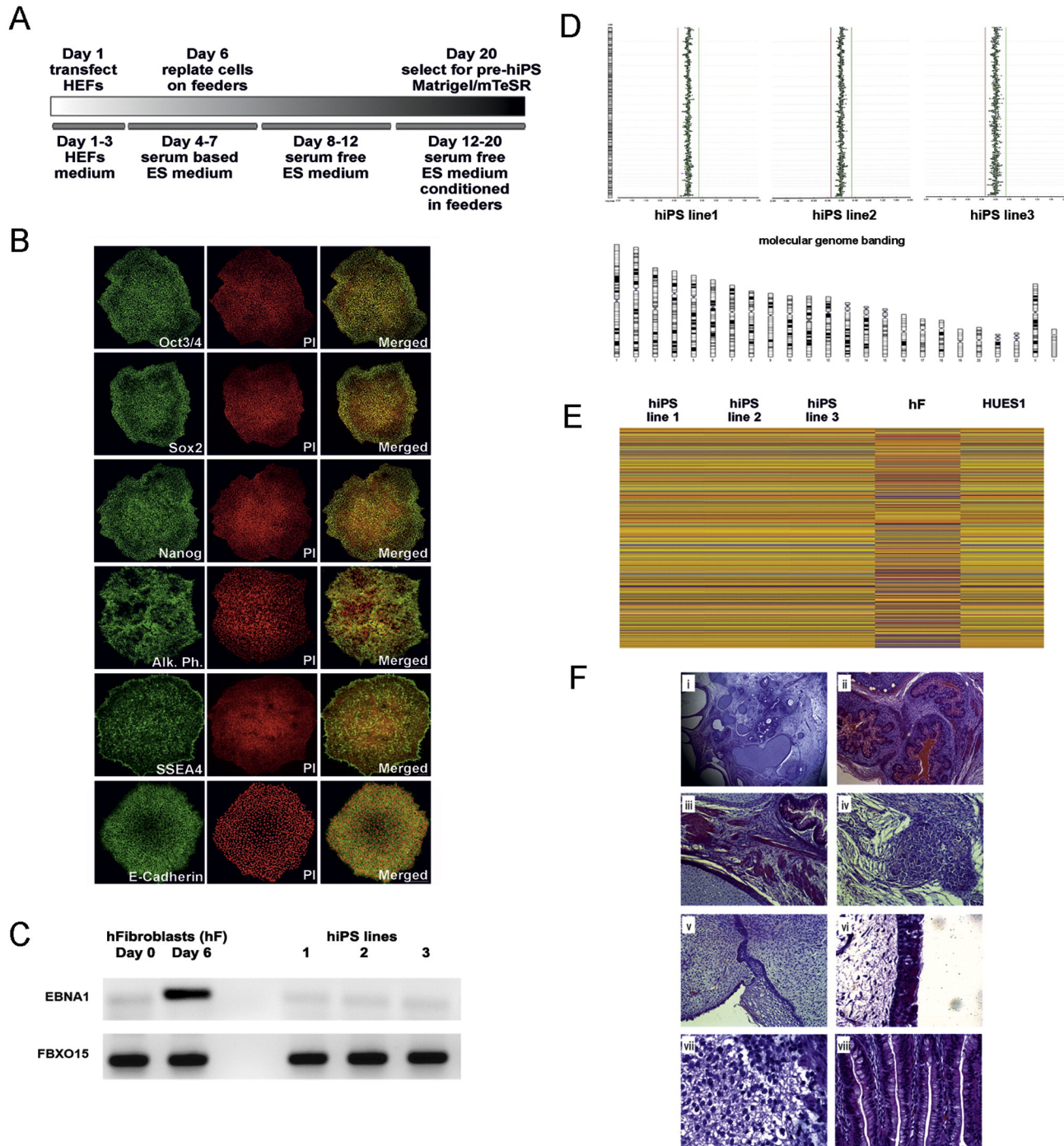


Fig. 1. Generation and validation of hiPS. A) Reprogramming strategy, B) immunofluorescence analysis of generated hiPS for embryonic stem cell markers. hiPS colonies were positive for Oct3/4, Sox2, Nanog, Alkaline phosphatase, SSEA4, and E-Cadherin stem cell markers. Propidium iodide was used to stain cell nucleus. C) Footprint analysis of generated hiPS lines. A PCR-based test for the validation of integration-free hiPS lines was performed targeting the EBNA-1 sequence derived from the Epstein-Barr virus to calculate the copy number of the episomal vector and another primer set for the endogenous locus FBXO15. The generated hiPS lines show no episomal footprint. D) Comparative genomic hybridisation analysis shows genetic stability with no evidence of chromosomal abnormalities, E) Heat map analysis of the differentially regulated genes following reprogramming among hiPSC lines, original human fibroblasts (hFs) and embryonic HUES1 stem cell line. All hiPS lines show heat map genetic profile identical to embryonic HUES and different from original hFs. F) Histological

markers (Fig. 1B). Cells were manually seeded (en bloc) on a matrigel coated chamber-slide system (LabTek #177437). Cells were fixed with either ice cold Methanol or paraformaldehyde 3.7% (followed by permeabilization with 0.1% Triton for 4 min). After blocking with 5% donkey serum, cells were stained with antibodies recognizing stem markers for 1 h at the following concentrations: Oct3/4 (1:300), Nanog (1:100), Sox2 (1:100), Alkaline Phosphatase (1:100), SSEA4 (1:100) and E cadherin (1:300) (Suppl. Table 1). After washing, secondary antibodies FITC donkey α -mouse or FITC donkey α -goat were added for 45 min, followed by 30 min incubation with 1 mg RNase (Sigma-Aldrich, Hertfordshire, UK) and 20 min with propidium iodide (Sigma-Aldrich, dilution 1:5000). Samples were mounted with Prolong antifade reagent (Invitrogen, Life Technologies) and imaging was performed with a Leica TCS SP5 (10 \times dry objective).

2.2. Genotyping

High-resolution array-based comparative genomic hybridization (CGH) (Centre of molecular biology and cytogenetics AlphaLab, Greece) showed that hiPSC lines had no indication for copy number variations and were therefore chromosomally stable (Fig. 1C). Moreover, PCR screening of the genomic DNA, using primers specific for the EBNA1 amplicon of the episomes used, verified footprint-free generation (Fig. 1D).

2.3. Whole transcriptome analysis

To validate the comparable transcriptional profile of the generated hiPSC lines versus human embryonic stem cells (HUES1 line, NIH repository), we proceeded with whole transcriptome analysis, accessible in the Gene Expression Omnibus (GEO) of NCBI

(#GSE58932). Quality-assessed RNA samples were hybridized to Affymetrix ExonExprChip.HuGene-1 arrays (EMBL-GeneCore, Heidelberg, Germany). Normalized data underwent statistical analysis (GeneSpring Software) to determine differentially expressed genes ($p \leq 0.05$ and fold change ≥ 2). A differential expression of approximately 3300 genes was observed between hiPSCs and control embryonic stem cell line HUES1 versus the initially transduced fibroblasts (heatmap in Fig. 1E).

2.4. Differentiation

To ensure full pluripotency we carried out teratoma assays in immune-deficient NOD/SCID mice, maintained under SPF conditions (Animal facility of the Academy of Athens, Greece). After manual dissection of hiPSC and HUES1 cultures (each from a subconfluent 35 mm dish), cell clumps were centrifuged and resuspended in 100 μ l of hES-qualified matrigel, just prior to injection. The cell-matrigel mix was injected through 25G 7/8 needle (BD Biosciences) subcutaneously in the hind legs of the mice. Transplanted animals developed tumors within a month and histopathological examination yielded tissue derivatives from all three germ lines: ectoderm, mesoderm and endoderm (Fig. 1F).

Supplementary data to this article can be found online at <http://dx.doi.org/10.1016/j.scr.2016.05.006>.

References

- Okita, K., Y. Matsumura, Y. Sato, A. Okada, A. Morizane, S. Okamoto, H. Hong, M. Nakagawa, K. Tanabe, K. Tezuka, T. Shibata, T. Kunisada, M. Takahashi, J. Takahashi, H. Saji and S. Yamanaka "A more efficient method to generate integration-free human iPS cells." *Nat. Methods* 8(5): 409–412.
- Kouroupis, D., Kyrkou, A., Triantafyllidi, E., Katsimpoulas, M., Chalepakis, G., Goussia, A., Georgoulis, A., Murphy, C., Fotsis, T., 2016. Generation of stem cell-based bioartificial anterior cruciate ligament (ACL) grafts for effective ACL rupture repair. *Stem Cell Res.* (submitted).